

# Luxor 250 Procedure Light

A sleek and high adjustable intensity LED procedure light with easy touch arm controls. This well-built line is designed to minimize contamination from blood, chemicals, and air circulation and is easy to clean. See color and detail like never before with an impressive 97 CRI.



## Illumination Intensity

3 Step adjustment (Max. 110,000 lx)



## Light Field Size

3 Step adjustment (Max. 25cm / Min. 15cm)



## Color Temperature

3,800K / 4,300K / 4,800K



## Color Rendering Index

Ra 97 / R9 97



ENDO

## Endo Mode

3 Step adjustment  
(Max. 15,000lx / Min. 5,000lx)



## Mounting Options

Single Ceiling  
Dual Ceiling  
Mobile  
Wall




**INFINIUM**<sup>®</sup>

SINCE 1991 

Exam . Oral . Dental . Derm  
ENT . Plastic . Cosmetic

# Bright, Reliable, & Intuitive


## TECHNICAL SPECIFICATIONS

Class	S250
Design (Light Head)	
Installation	Single Ceiling Dual Ceiling Mobile Wall
Head Size (cm)	37.2 X 38
No. of LED (ea)	34
Maximum Intensity (lx)	110,000
Intensity control	3 Step (40,000 ~ 110,000)
C.R.I (Ra / R9)	97 / 97
Color Temperature (K)	3,800 / 4,300 / 4,800
Light Field Size (Cm)	3 Step (15 ~ 25, Circle)
Detachable Centre grip	Yes (2X)
Input Power	100 ~ 240Vac, 50/60Hz
LED Life	60,000Hrs
Option	Backup Battery


## Luxor 250 Surgical Lights

Premium procedure lights are priced to fit your budget. The Luxor 250 series is long-lasting, feature-rich, easy to use, and backed by a 3yr warranty

### Easy Touch Control Systems

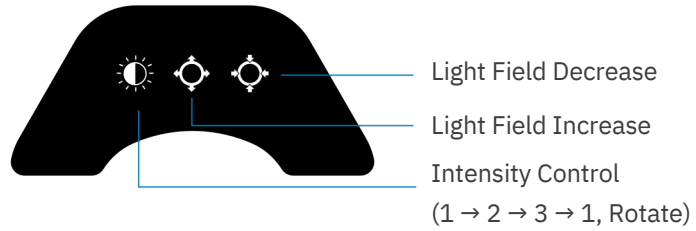


### ARM Touch Control (Intuitive)

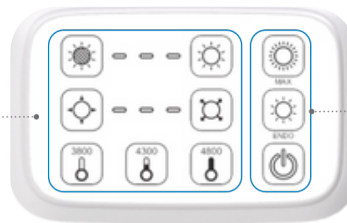

### Light Control

- Brightness (Intensity) 3 Step
- Light Field Size 3 Step
- Color Temp. 3 Step

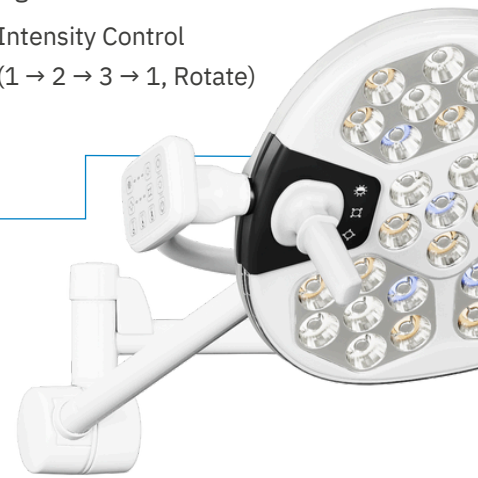

### Max. Mode

- Light Field Size 3 Level -
- Intensity 3 Level

### Endo Mode

- Light Field Size 1 Level -
- Intensity & Color Temp. 3 Step (5,000 ~ 15,000 lx)

### Light On / Off



# INFINIUM®

Infinium Medical  
12151 62nd St North #5  
Largo, FL 33773 • USA

(1) 866-918-8434  
infiniummedical.com  
sale@infiniummedical.com