



ATS TRINITY

CARBON FIBER TOP IMAGING TABLE

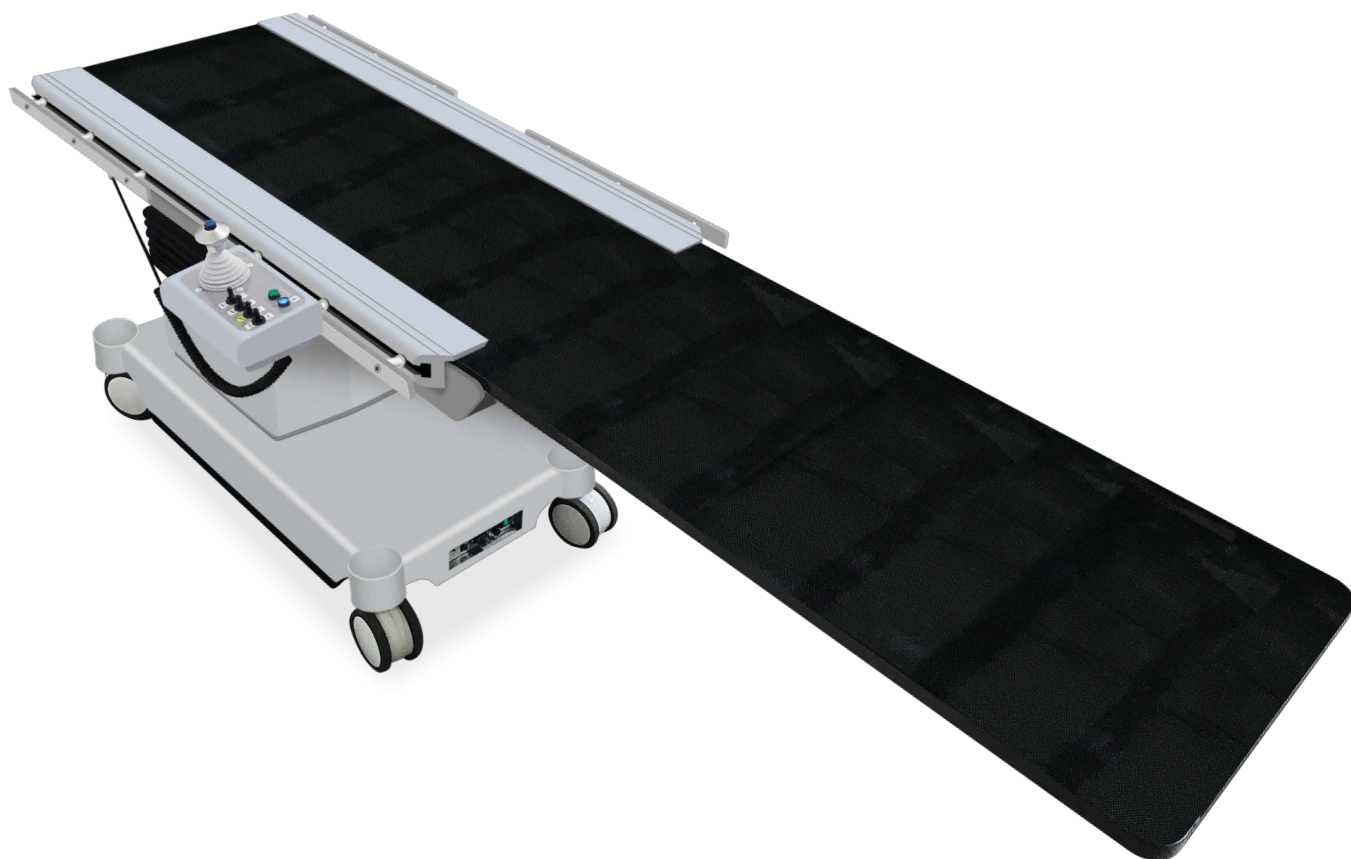

INTUITIVE MOTION
UNYIELDING STABILITY

ATS LATERIS

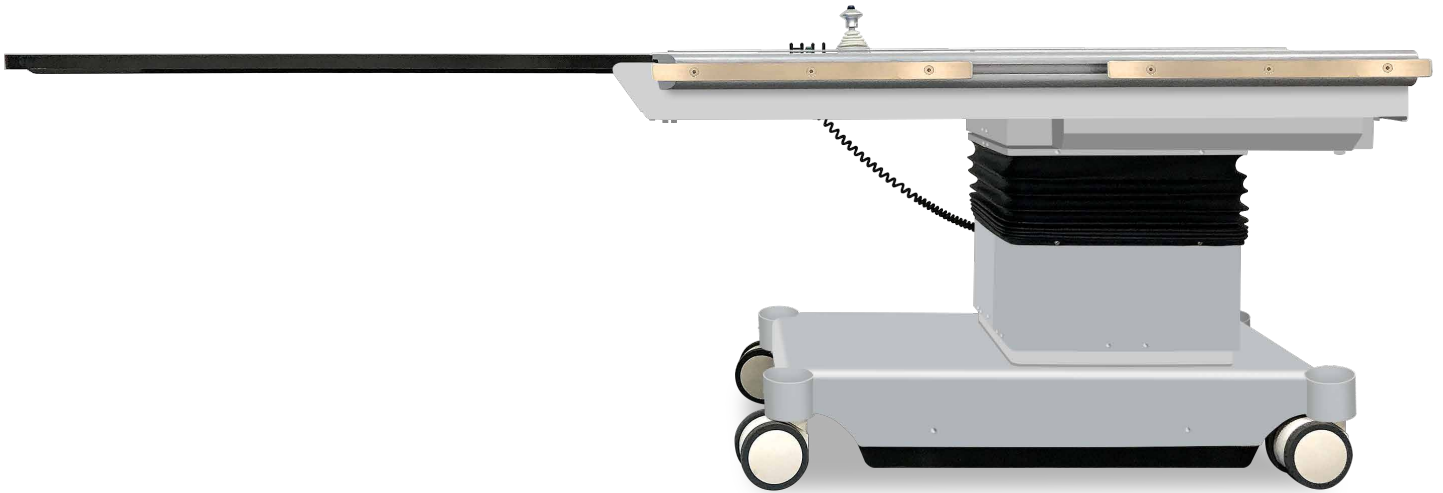

Intuitive Motion

Built to the foremost standard of quality, with intuition and stability in mind. The ATS Trinity® imaging table exceeds the standards in patient positioning. The ATS surgical table features a C-arm compatible carbon fiber table top with a 550 LB lift capacity. The ATS is equipped with a motorized 4 way floating top to maximize patient imaging access. Other movements include: Trendelenburg, Reverse Trendelenburg, Lateral Tilt, and Up/Down elevation.

Patient Safety


The ATS Surgical table is mobile on 4 anti-static castors. A floor lock system braces the table in place with unyielding stability. A battery backup system is a standard feature. The removable hand controller is ergonomic and outfitted with safety locks to avoid unintentional table movement. A return to level function is available as an option.

ATS Trinity Features

- Carbon Fiber Top
- Silent Electro-Hydraulic Movements
- 550 LB Lift/Weight Capacity
- Electronic Control Safety Lock
- Mobile Anti-static Castors with Floor Lock
- Floating Top
- Left and Right 15 Degree Lateral Lilt
- Trendelenburg / Reverse Trendelenburg
- Battery Backup
- Accessory Rails
- Optional Auto-Reset / Level Function

Built To The Foremost Standard Of Quality


WITH MOVEMENT AND STABILITY IN MIND.

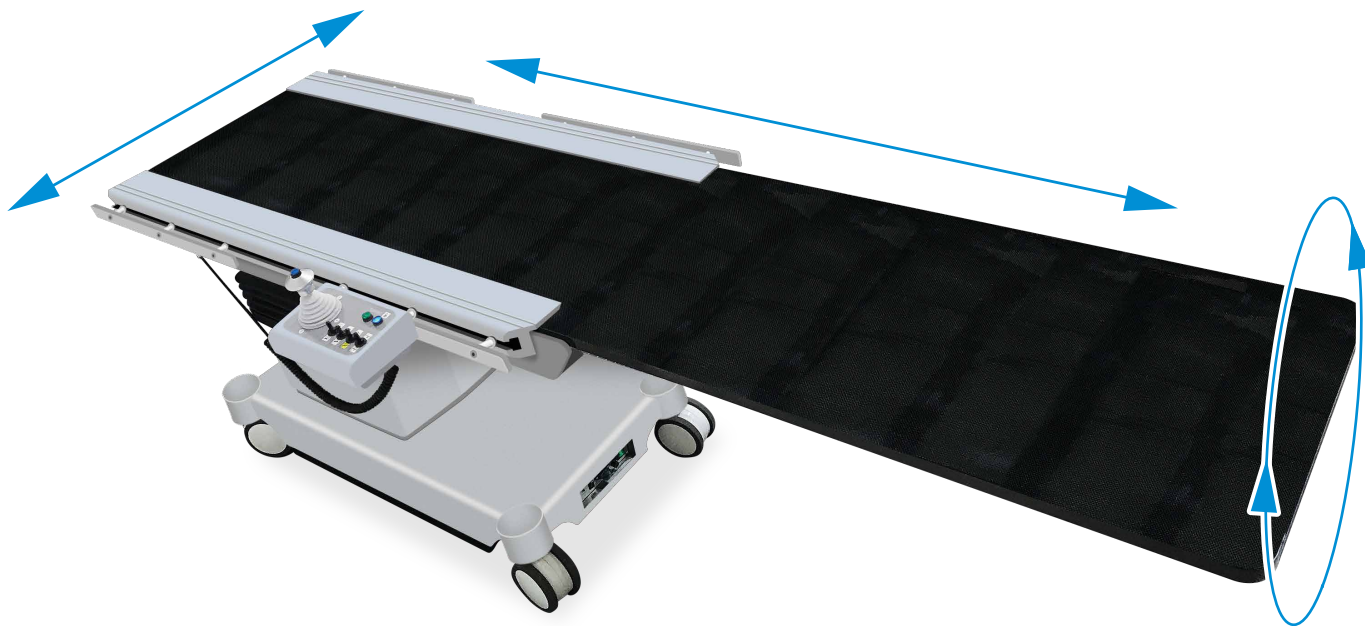

- 4 Way Floating Top with Lateral Tilt




- Intuitive All-In-One Removable Hand Control

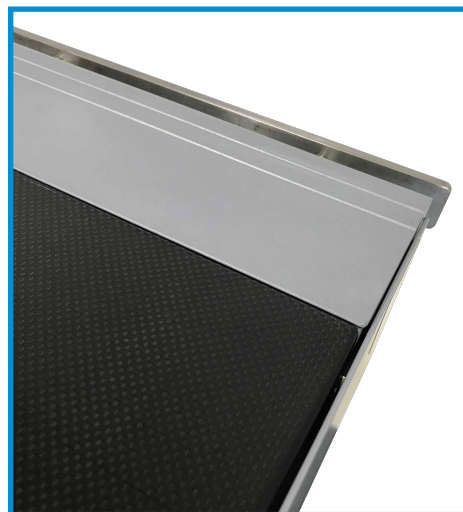


- Carbon Fiber Top with Removable Accessory Rail Mounts



ATS TRINITY TECHNICAL SPECIFICATIONS:

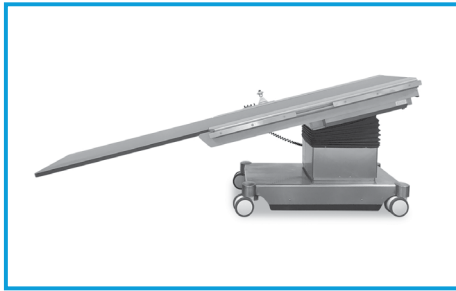



■ Trendelenburg

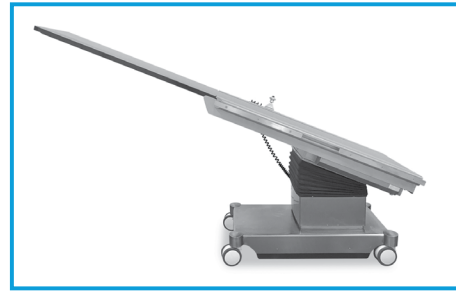

■ Reverse Trendelenburg

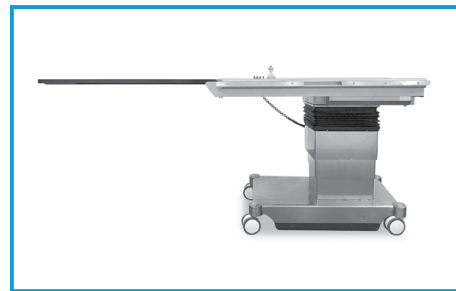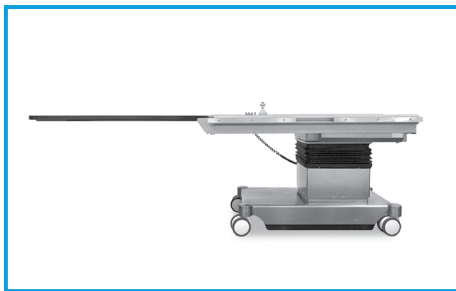

■ Adjustable Elevation

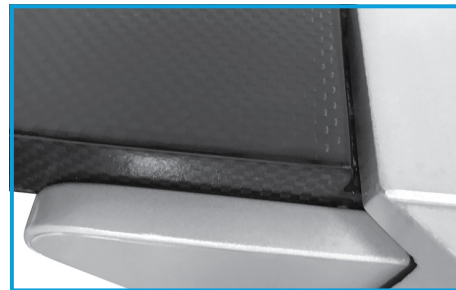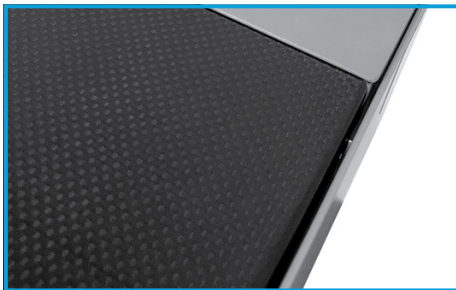

■ Carbon Fiber and Stainless Steel Construction

SPECIFICATIONS:

Table Top Length	2500 mm
Table Top Width	500 mm

MOVEMENTS:

Adjustable height	700 mm - 1100 mm
Trendelenburg / Reverse Trendelenburg	30 Degrees
Oblique left/Right Lateral Tilt	15 Degrees

Longitude Table Top Float	700 mm
Side/Side Table Top Float	200 mm
Weight Capacity	250Kg (550 LBS)
Ingress Protection	IPX4
Power Supply	Universal AC 110V-220V / 50Hz-60Hz 50Hz-60Hz DC 24V
Battery Backup	2.4 AH